
# A E N O


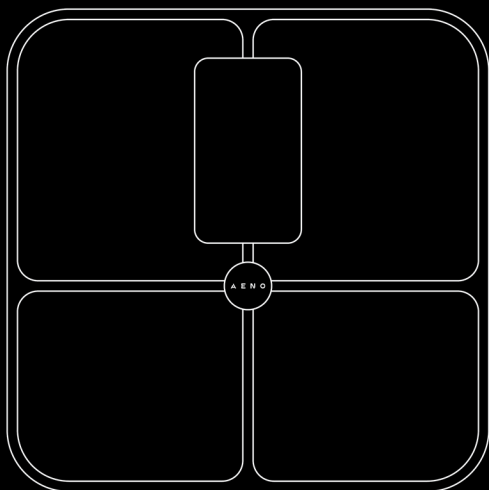
Quick Start Guide

## Smart Body Scales BS2S

(Model: ABS0002S)

ENG	4	KAZ	18
ARA	5	LAV	19
BOS	6	LIT	20
BUL	7	NLD	21
CES	8	POL	22
DEU	9	POR	23
ELL	10	RON	24
EST	11	RUS	25
FRA	12	SLK	26
HRV	13	SLV	27
HUN	14	SPA	28
HYE	15	SRP	29
ITA	16	UKR	30
KAT	17	UZB	31

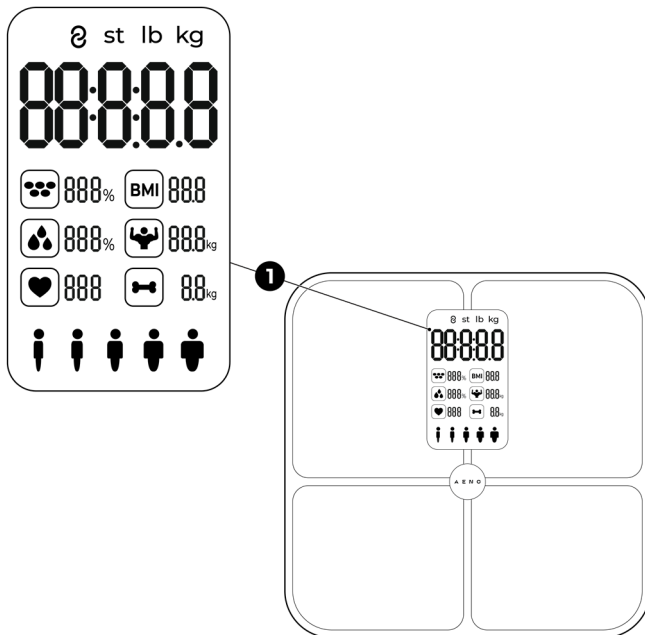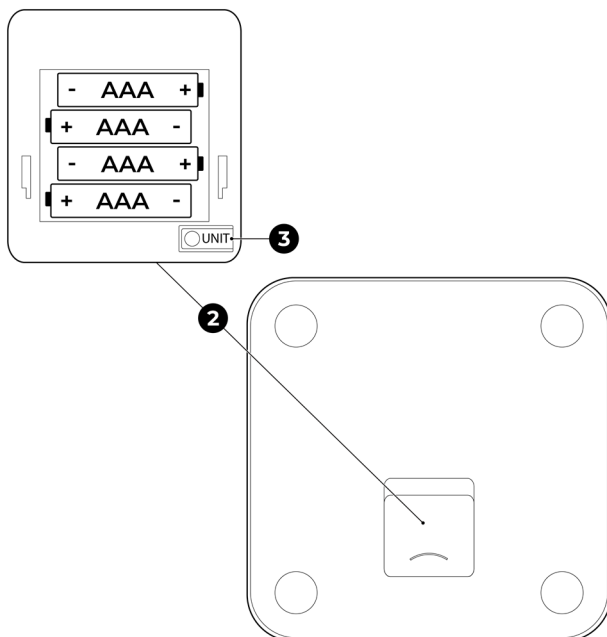**A****B**

**ENGAENO™ ABS0002S** Smart Body Scales are designed to determine the user's weight as well as additional body parameters.

**Technical Specifications**

Power supply: AAA batteries (4 pcs.), 6.0 V, 0.4 W (max.). BT 5.0. Operating frequency: 2,400–2,483 MHz. Transmitter power: +5 dBm / 3.16 mW. Measuring range: 5–180 kg. Weight increment: 0.1 kg. Measurement error: ±0.3 kg (for 0–50 kg); ±0.43 kg (for 50–100 kg); ±0.5 kg (for more than 100 kg). Units of measure: kg, lb, st/lb. Display: VA; 107×62 mm. Automatic activation/deactivation. Indication of the BT connection to the app. Calculation of body fat percentage, muscle and bone mass, body water level, heart rate, body mass index, body type. Control: manual, remote (AENO app). Color: Black. Casing material: HIPS plastic, tempered glass with ITO coating. Dimensions (L×W×H): 278×278×25.5 mm. Weight: 1.4 kg. Operating conditions: temperature +5...+40 °C, RH 40–85 % (non-condensing). Storage conditions: temperature –20...+60 °C, RH up to 90 % (non-condensing).

**Scope of Supply**

Body Scales **ABS0002S**, AAA batteries (4 pcs.), Quick Start Guide.

**Restrictions and Warnings**

Handle the glass surface of the scales with care: avoid damaging, cracking and chipping the surface as broken glass can cause injury. Supervise the use of the device by children and persons with physical and mental disabilities. Do not allow children to play with the device. Place the scales on a hard, flat, horizontal surface to avoid inaccurate measurements. Do not place the scales on carpets or any other soft surface. Do not load the device beyond the maximum value of the measuring range. If the device is not used for an extended period of time, remove the batteries and store them out of the reach of children and animals. Do not use chemical and harsh detergents, abrasive pastes, products containing acids and solvents, or metal sponges to clean the device.

**ATTENTION!** Do not expose the device to water or other liquids. Do not immerse the device in water. Do not wash the device or any part of it under running water or in a dishwasher.

**Device Elements** (see Fig. A, B)

1 – display; 2 – battery compartment; 3 – “UNIT” button.

**Device Indicators**

	Active pairing with the AENO app		Heart rate
	Selected weight measurement unit		Body mass index value
	Weight display scale		Muscle mass level
	Body fat percentage		Bone mass level
	Body water level		Body type of the user

**Preparing the Device for Operation**

Unpack the device and clean it (see “Cleaning and Maintenance”). Open the battery compartment (2) cover and insert the batteries into the compartment, observing the correct polarity. Press the “UNIT” button (3) to select the unit of weight measurement. After each press, the display (1) will show the current unit of measurement: “kg” – kilogram, “lb” – pound, “st/lb” – stone. Close the cover of the compartment. Place the device on a hard, horizontal surface.

**Note.** The selected settings are saved for current use and for the next time the device is turned on. To reset the settings, press and hold the “UNIT” button for 5 seconds. The display will show the code “CLr”. The settings are also automatically reset when the batteries are removed.

**Operation of the Device**

Stand on the surface of the device. The scales will turn on automatically, and the display will show your weight in the selected unit of measure. Within 5–10 seconds after the weight measurement is complete, the display will show the remaining measurements. Note that the user must stay on the surface of the scales. The scale will turn off automatically within 15 seconds after all measurements have been completed.

**ATTENTION!** Information shown by the scales cannot be considered medical information.

**Cleaning and Maintenance**

**Note.** Remove the batteries before cleaning the device.

Wipe the device with a soft cloth dampened with clean water and then wipe dry. If the scales are very dirty, wipe them with a soft cloth dampened with a mild soap solution. Then wipe with a soft cloth dampened with clean water and wipe dry.

**Connecting to AENO Mobile App**

Download the AENO mobile app from Google Play or App Store and register a new account (or log in to an existing one).

- Stand on the scales to turn them on. Place your smartphone as close to the device as possible.
- On the home screen of the app, tap “+” and then tap “Add Device”.

**Note.** If you receive system prompts to turn on any features on your smartphone, follow the prompts.

- The device will be found automatically, and its image will appear in the “Autosearch” section. Tap the “Add” button.
- Wait until the message indicating that the device has been successfully added appears. Tap “Done”.

All parameters that are shown on the display of the scales are available in the app, as well as information about the user's metabolic rate and biological age, physical condition rating, etc. For more information on the features available in the app, refer to the complete user manual at [aeno.com/documents](http://aeno.com/documents).

**Troubleshooting**

- The device does not turn on. Possible cause: batteries are missing or discharged; polarity was not observed when inserting the batteries. Solution: insert charged batteries into the compartment, observing the polarity.
- The device does not connect to the app. Possible cause: BT is not enabled on the smartphone; the smartphone is too far away from the scales; scales settings failure. Solution: enable BT on your smartphone; place the smartphone no more than 5 m away from the scales; reset the scales to factory settings.

**Error Codes Shown on the Display**

“Lo” – low power supply charge level; “Err” – maximum allowable weight exceeded; “C” – weighing failure.

**ATTENTION!** If none of the suggested methods have solved the problem, contact your supplier or an authorized service centre. Do not disassemble the device or attempt to repair it yourself.

If you have any questions or difficulties in using your AENO™ device, please e-mail support at [support@aeno.com](mailto:support@aeno.com) or chat with us at [aeno.com/service-and-warranty](http://aeno.com/service-and-warranty). The specialists will help you solve the problem, and you will not have to waste time and effort visiting the store.

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**ENG** ASBISIC reserves the right to modify the device and to make changes and additions to this document without prior notice to users.

The warranty period and service life is 2 years from the date of retail sale of the device.

Manufacturer: ASBISIC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Made in China. All trademarks and brands mentioned herein are the property of their respective owners. The date of manufacture is indicated on the package.

Up-to-date information and a detailed description of the device, as well as connection instructions, certificates, information on companies that accept quality claims and warranties are available at [aeno.com/documents](http://aeno.com/documents).

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**ENG Recycling information.** These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility. You should contact your local household waste disposal service for details.

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