A E N O

Quick Start Guide

Dental Irrigator


D11

(Model: ADI0001)


10.10.2023 | v2.0.0

ENG.....	13
ARA.....	17
BUL.....	20
CES.....	24
DEU.....	28
ELL.....	32
EST.....	36
FRA.....	40
HRV.....	44
HUN.....	48
HYE.....	52
ITA.....	56
KAT.....	60
KAZ.....	64
LAV.....	68
LIT.....	72
NLD.....	76
POL.....	80
POR.....	84
RON.....	88
RUS.....	92
SLK.....	96
SPA.....	100
SRP.....	104
UKR.....	108
UZB.....	112
ZHO.....	116


A

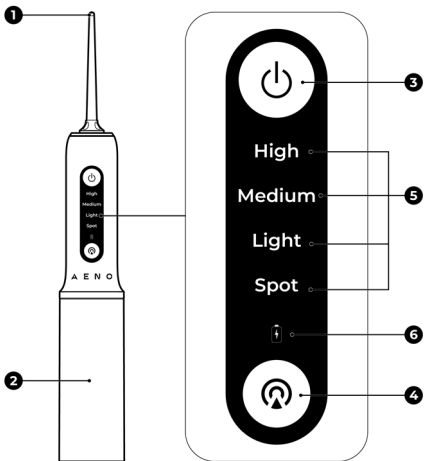
B

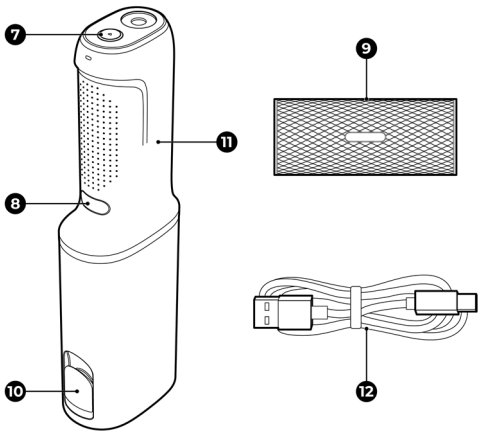
C

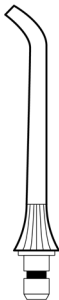13

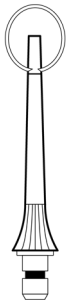14

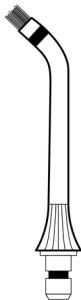
15



16

D


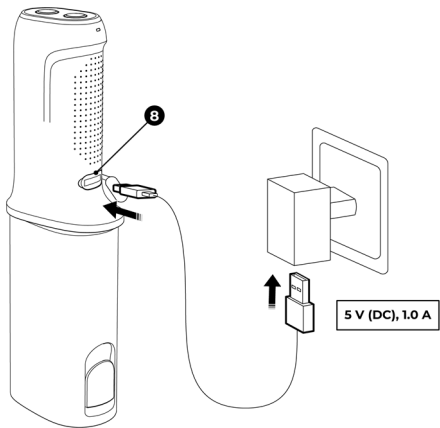
E

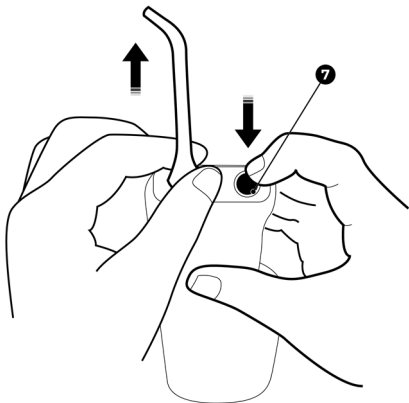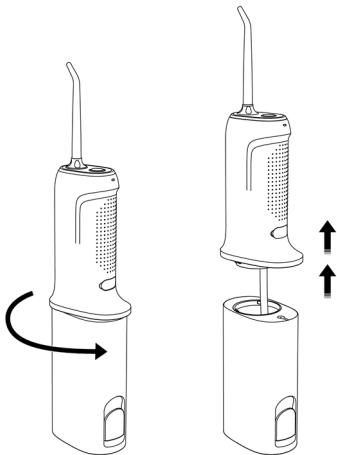
F

G

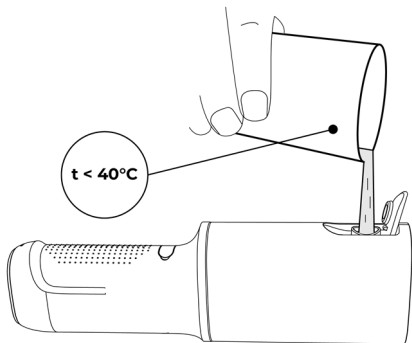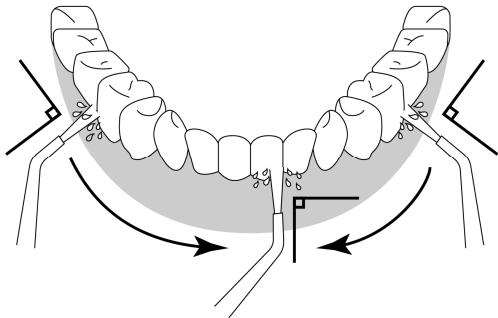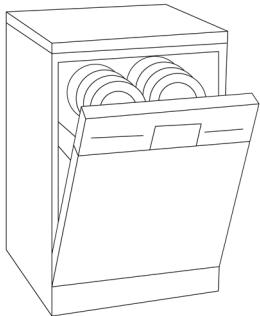ENG AENO™ ADI0001 Irrigator is designed for oral care, cleaning plaque and food debris in hard-to-reach places.

Specifications

Power input parameters: 5 V (DC), 1.0 A, 5.0 W (max.) Battery: non-detachable, lithium-ion; 2,200 mA-h, 3.7 V (DC). Noise level (at distance of 0.3 m): ≤ 69 dB. Water pressure: 206–689 kPa. Pulsation: 1,500 pulses/min. Full battery charge time: 4.5 h. Water tank: 200 ml. Casing protection degree: IPX7. Operating time after a full charge: no less than 45 minutes. Operating time with a full tank: no less than 55 s. Dimensions (L×W×H): 70×45×223 mm. Operating conditions: temperature +5...+40 °C, relative humidity 45–75 %. Storage conditions: temperature 0...+50 °C, relative humidity 45–90 %.

Package contents (figures B, C)

AENO™ ADI0001 irrigator (white) (11), 4 replaceable nozzles (13–16), USB (A)/USB (C) cable (1 m) (12), case for replaceable nozzles (9), Quick Start Guide.

Limitations and Warnings

Monitor the use of the device by children under the age of 12 and persons with physical or mental disabilities. Do not allow children to play with the device. If you have oral diseases, you should consult your dentist before using the device. Do not turn on a damaged device. Do not leave the plugged-in device unsupervised. Do not use the irrigator and do not perform maintenance on the device or its accessories unless the device is disconnected from the power supply. **Please use purified water without adding mouthwash, or other liquids or substances.** Do not point the nozzle outlet into your eyes while the device is operating.

If you have any questions or issues with your AENO device, please contact our online support team support@aeno.com or aeno.com/service-and-warranty. This is a better option than visiting the store and will save you time and effort.

Preparing the device for operation

Unpack the device. Wipe its casing with a cloth soaked in a little water and wipe it dry. Disconnect the tank (see fig. **F**). Rinse and dry the tank and nozzles as described in "Maintenance" section, clause "Daily Care".

Charging the device

WARNING! The power adapter is not included in the package. You can use the AC adapter from a smartphone, tablet or other device as long as it provides 5 V (DC), 1.0 A power parameters.

Connect the device to the mains as shown in figure **D**. Status of the battery charge indicator (6) during charging:

- the device is charging – the indicator blinks red;
- the device is fully charged – the indicator lights green.

When charging is completed, unplug the device, disconnect the cable (12) and close the USB port (8) by firmly inserting the plug. Charge the device for at least 12 hours before using it for the first time.

Status of the battery charge indicator (6) in the operating mode:

- the charge level is less than 5 % – the indicator blinks red rapidly, the device does not turn on;
- the charge level is 5–10 % – the indicator blinks red slowly;
- the charge level is 10–100 % – the indicator is off.

When the device is charged, place one of the replaceable nozzles (13–16) on the casing by pushing it until it clicks. The nozzle can spin 360°.

Purpose of the irrigator nozzles (figure **C**):

- standard nozzle for cleaning food debris, plaque (13);
- tongue cleaner nozzle (14);
- orthodontic nozzle for braces, dentures (15);
- nozzle for cleaning of the periodontal pockets (16).

Fill the water tank (2) as shown in figure **G**. Close the water fill cap (10) tightly. After that, go on to select the operating mode of the device.

Setting the operating mode of the device

To set the operating mode, point the nozzle's outlet in a safe direction so that the water jet does not damage the interior and so that you can see the mode selection panel (5).

Press the button (3) to turn on the device. The operating mode indicator on the irrigator panel (5) will light up. Use the mode selection button (4) to select the desired mode. The first time you turn on the device, the default operating mode is "High". Subsequent single presses of the mode selection button (4) alternately change the operating modes, which are displayed on the panel (5).

Device operating modes

- "High" – intensive cleaning for healthy gums;
- "Medium" – daily standard cleaning, plaque removal;
- "Light" – gentle cleaning for sensitive teeth and gums. Also recommended for first-time users of the device;
- "Spot" – removal of food debris from interdental spaces, prevention of tartar formation. Water is ejected when button (3) is pressed.

Each mode of the device has three power levels: high, medium, low. To switch the power level, press and hold the button (4) for 3 (three) seconds. The mode and its power level, set during the last use, will be automatically selected the next time the irrigator is turned on. Turn off the device by using the button (3).

Device operation

Point the nozzle outlet (1) into your mouth cavity. Press the button (3) to turn on the device. Move the nozzle's outlet (1) along the gum contour as shown in figure **H**. During the cleaning, you can switch the operating mode and its power level with the button (4). When you have finished the cleaning, press the button (3) to turn the device off. If you do not turn off the device, it will automatically turn off not later than in 2 (two) minutes. After each use, wash the device and the used nozzles as indicated in the "Maintenance" section (clause "Daily Care").

NOTE. It is not recommended to use the nozzle (16) in the "Light" mode to avoid damage to periodontal pockets.

Maintenance

- Daily Care. Open the filler cap (10) and drain off any water remaining in the tank. Disconnect the nozzle as shown in figure **E**. Rinse the tank and nozzles in warm running water, without detergents or cleaning agents, and dry thoroughly.
- Total care. A total cleaning of the water tank must be performed once a week. To do this, disconnect the tank as shown in picture **F**. Wash it inside and outside in warm running water with a soft brush or sponge, without detergents or cleaning agents, and dry thoroughly.

Troubleshooting

- The device does not turn on. Possible reason: the battery is low. Solution: charge the device for at least 4.5 hours.
- The device does not supply water. Possible causes: the nozzle's outlet is clogged, the water tank is out of water. Solution: fill the tank with water, clean the nozzle outlet.
- Operating time after a full charge is less than 10 minutes. Possible reason: the battery does not hold charge. Solution: contact the service centre to replace the battery of the device.

WARNING! If none of the suggested solutions solves the problem, contact your supplier or service centre. Do not disassemble or attempt to repair the device yourself.

ASBISc reserves the right to modify the device and to make changes and additions to this document without prior notice to users.

The warranty period and service life is 2 years from the date of retail sale of the product.

Manufacturer's information: ASBISc Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Made in China. All trademarks and names herein are the property of their respective owners. Up-to-date information and a detailed description of the device, as well as connection instructions, certificates, information about companies that accept quality claims and warranties, are available for download at aeno.com/documents.

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ENG Recycling information These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to

a local recycling facility. You should contact your local household waste disposal service for details.

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