
MOVA Electric Toothbrush

User Manual

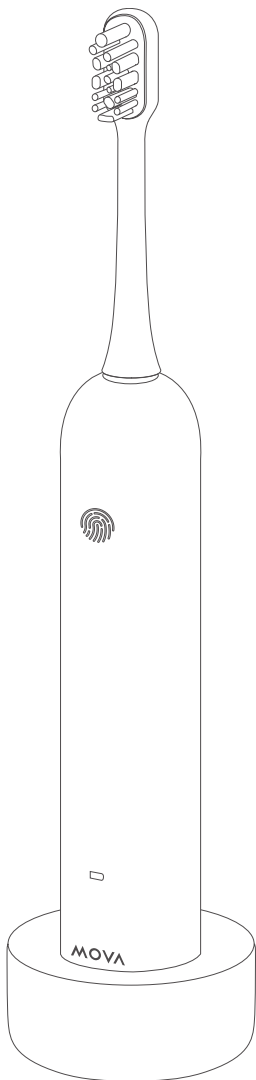*The illustrations in this manual are for reference only.
The actual appearance of the product may vary from
what is shown.

MOVA

EN	User Manual.....	001
DE	Benutzerhandbuch.....	006
FR	Manuel d'utilisation	011
IT	Manuale dell'utente	016
ES	Manual de usuario	021
NL	Gebruikershandleiding.....	026
NO	Brukerhåndbok	031
SV	Användarhandbok	036
FI	Käyttöohje.....	041
DK	Brugermanual.....	046
PL	Instrukcja obsługi.....	051
TR	Kullanıcı Kılavuzu	056
TH	คู่มือการใช้งาน.....	061
VN	Hướng Dẫn Sử Dụng.....	066
ID	Panduan Pengguna	071
HE	מדריך למשתמש.....	076
AR	دليل المستخدم	081
KO	사용 설명서	086
TW	使用說明書	091
JP	取扱説明書	096
PT	Manual do utilizador	101
GR	Εγχειρίδιο χρήστη	106
KZ	Пайдаланушы нұсқаулығы	111
UA	Посібник користувача	116
GE	მომხმარებლის სახელმძღვანელო	121
UZ	Foydalanuvchi qo'llanmasi	126

Scan the QR code to download the MOVAhome APP and connect to the electric toothbrush for more functions. Before connecting the electric toothbrush to your phone for the first time, please turn on the electric toothbrush at least once to ensure it is activated.

Scan the QR code to view the product tutorial and electronic manual

Product Description

Packing List

Electric toothbrush handle	x1
3-in-1 brush head	x1
Sensitive care brush head	x1
Wireless charging base	x1
User manual	x1
Safety tips	x1
Travel case	x1

Part Description

Brush head

3-in-1 brush head

Triple functions including gentle tooth protection, efficient cleaning, whitening and polishing; suitable for all people

Sensitive Care Brush Head

Ultra-fine and ultra-soft brush bristles are adopted, reaching deep into tooth gaps and caring for gums, suitable for people with sensitive teeth

Status indicator

Display different colors in different brushing modes

Power/Mode switching button

Battery/Bluetooth indicator



Dust-proof cover to protect brush head

Multi-functional travel box for brush head

Used to store the brush head when traveling

Travel case

Wireless charging base

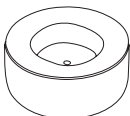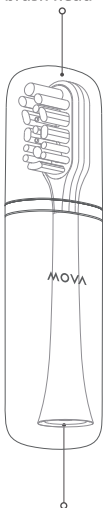
Instructions for Use


EN

- ① Insert the brush head into the handle
- ② First Use: Long press to activate
When the toothbrush is not working, long press the button to enter/exit Flight Mode.

* Toothbrush beeps 3 times when entering or exiting Flight mode.

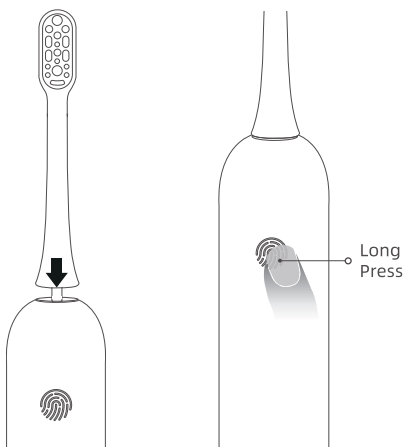- ③ Short press to turn on/ switch modes.
Long press to pause.

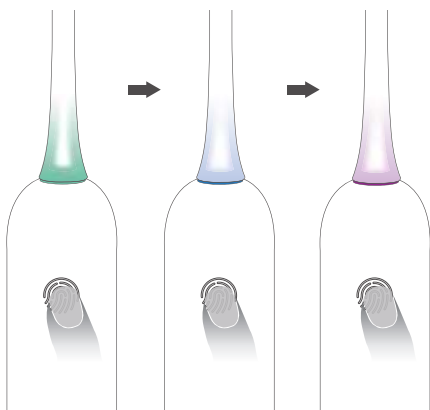Green - Sensitive Gentle Mode
Suitable for first-time use people with sensitive teeth


Blue - Standard Cleaning Mode
Suitable for daily use and people seeking efficient cleaning

Purple - Whitening & Polishing Mode
Suitable for deep cleaning people who want brighter teeth

Smart Functions

Smart pressure control: If excessive brushing pressure is detected during use, the status indicator will flash in red and the toothbrush will automatically reduce its vibration frequency to prevent potential tooth damage.

Flight Mode: Entering Flight Mode can prevent toothbrush from being activated by short accidental touches while traveling. In Flight Mode, the light ring glows white when the button is pressed. When the toothbrush is not working, long press the button to enter or exit Flight Mode.

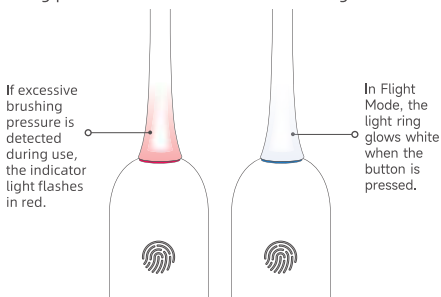Tip: The above functions can be enabled or disabled in the MOVAhome APP.


Battery

Battery Level Reminder

Battery Level	40%-100%	20%-39%	10%-19%	<10%
Battery Indicator	Solid green	Solid orange	Solid red	Flashing in red
				
	MOVA	MOVA	MOVA	MOVA


Charging

1. Connect the power cord of the charging base to the power adapter.
2. Place the toothbrush handle correctly onto the charging base.
3. Do not immerse the charging base in water or any liquid.
4. Do not connect the power source to the charging base if there is any liquid on the power plug.

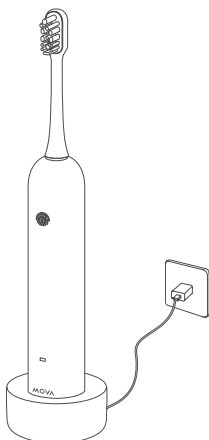Indicator status:

Battery level 0-99%	: Breathing in green
Battery level 100%	: Solid green
During Bluetooth connection	: Breathing in blue
Bluetooth is successfully connected	: Solid Blue

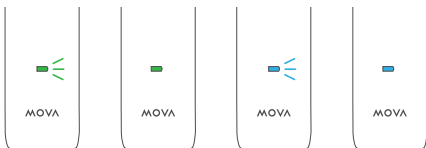Basic Parameters

Product name	Electric toothbrush		
Model	ATB15A		
Rated input voltage	5V 	Rated input power	5W
Waterproof rating	IPX7	Battery capacity	900mAh

The Bluetooth LE:

max e.i.r.p. : 2.3dBm

Frequency band of operation: 2400-2483.5MHz

Wireless charger :

Frequency Range: 110-205kHz

H-field : 60.1dBµV/m @3m

Operating Temperature: -1

0 °C ~ 55 °C

Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.

Children shall not play with the appliance.

This appliance contains batteries that are non-replaceable. Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmentally safe recycling.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

This appliance contains batteries that are non-replaceable.

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Manufacturer: Moonstone Technology (Shaoxing) Co., Ltd.