MOVA Electric Toothbrush

User Manual

*The illustrations in this manual are for reference only.
The actual appearance of the product may vary from
what is shown.

MOVA

EN	User Manual.....	001
DE	Benutzerhandbuch.....	006
FR	Manuel d'utilisation	011
IT	Manuale dell'utente	016
ES	Manual de usuario	021
NL	Gebruikershandleiding.....	026
RU	Руководство пользователя.....	031
NO	Brukerhåndbok	036
SV	Användarhandbok	041
FI	Käyttöohje.....	046
DK	Brugermanual.....	051
PL	Instrukcja obsługi.....	056
TR	Kullanıcı Kılavuzu	061
TH	คู่มือการใช้งาน.....	066
VN	Hướng Dẫn Sử Dụng.....	071
ID	Panduan Pengguna	076
HE	מדריך למשתמש.....	081
AR	دليل المستخدم	086
KO	사용 설명서	091
TW	使用說明書	096
JP	取扱説明書	101
PT	Manual do utilizador	106
GR	Εγχειρίδιο χρήστη.....	111
KZ	Пайдаланушы нұсқаулығы	116
UA	Посібник користувача	121
GE	მომხმარებლის სახელმძღვანელო	126
UZ	Foydalanuvchi qo'llanmasi	131

Scan the QR code to view
the product tutorial and
electronic manual

Product Description

Packing List

Electric toothbrush handle	x1
3-in-1 brush head	x1
Sensitive care brush head	x1
Charging cable	x1
User manual	x1
Safety tips	x1

Part Description

Brush head

3-in-1 brush head

Triple functions including gentle tooth protection, efficient cleaning, whitening and polishing; suitable for all people

Sensitive Care Brush Head

Ultra-fine and ultra-soft brush bristles are adopted, reaching deep into tooth gaps and caring for gums, suitable for people with sensitive teeth

Status indicator

Display different colors in different brushing modes

Power/Mode switching button

Battery/Bluetooth indicator

Dust-proof cover to protect brush head

Multi-functional travel box for brush head

Used to store the brush head when traveling



Charging cable

Instructions for Use


① Insert the brush head into the handle

② Short press the button once to turn on/pause

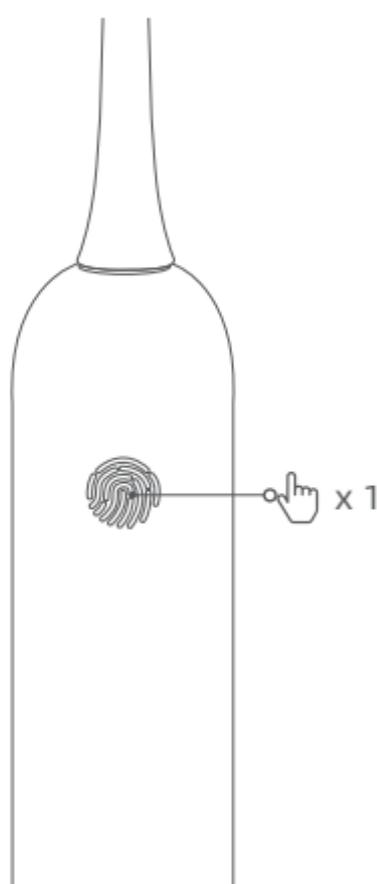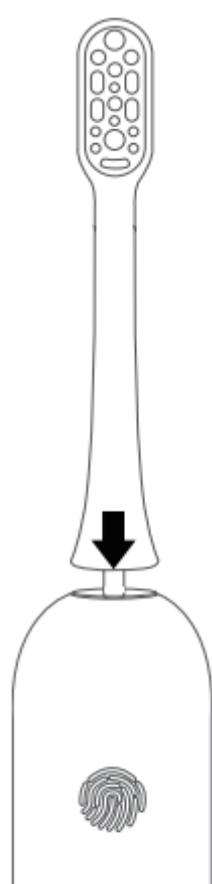Mode Switching

After turning on, press the button to switch modes within the first 8 seconds, and press again to pause after it has been working for 8 seconds.

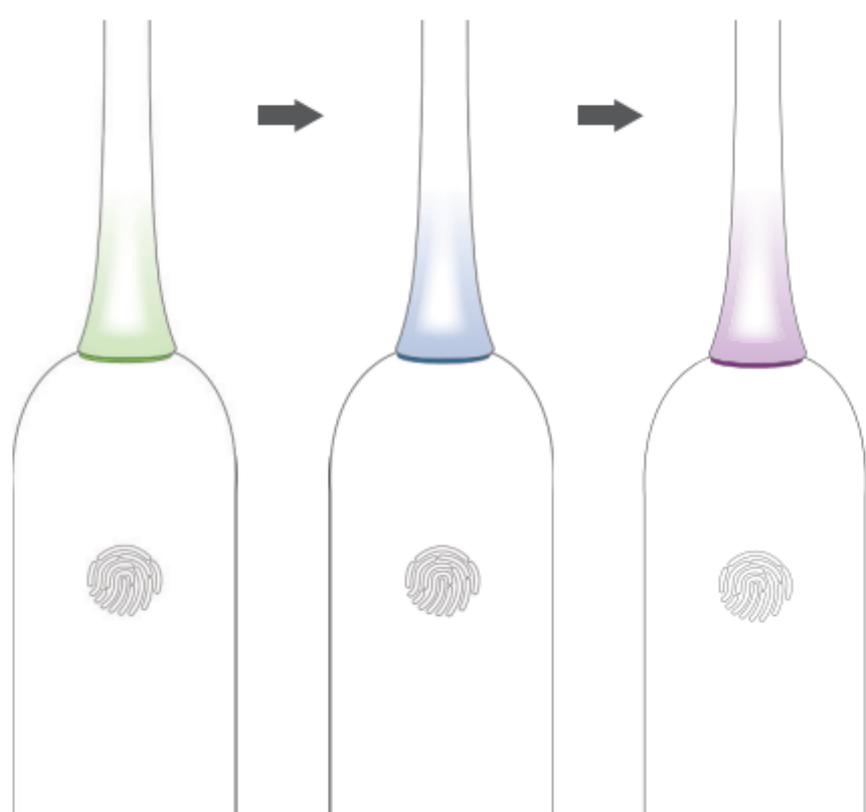Green - Sensitive Gentle Mode
Suitable for first-time use People with sensitive teeth

Blue - Standard Cleaning Mode
Suitable for daily use and people seeking efficient cleaning


Purple - Whitening & Polishing Mode
Suitable for deep cleaning People who want brighter teeth

Smart Functions

Smart pressure control: If excessive brushing pressure is detected during use, the status indicator will flash in red and the toothbrush will automatically reduce its vibration frequency to prevent potential tooth damage.

Smart mistouching protection: The toothbrush can only be activated when the button is directly pressed by the user's finger. The toothbrush is designed to prevent accidental activation during transportation, even if other items press against the button.

If excessive brushing pressure is detected during use, the indicator light flashes in red.


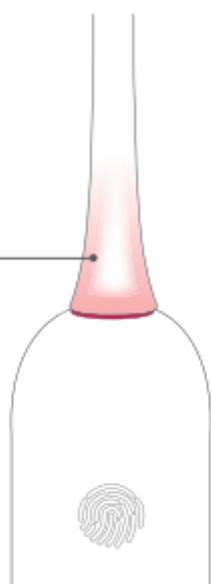
Battery

Battery Level Reminder

Battery Level	40%-100%	20%-39%	10%-19%	<10%
Battery Indicator	Solid green	Solid orange	Solid red	Flashing in red
				
	MOVA	MOVA	MOVA	MOVA

Charging


1. Open the charging cover at the bottom of the toothbrush and insert the charging cable correctly.
2. Use a USB charging adapter that complies with national standards.

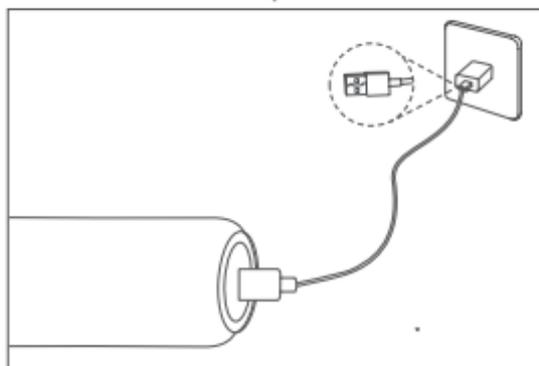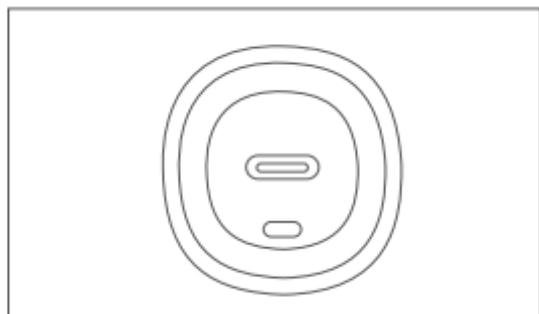Indicator status during charging:

Battery level 0-99%: Breathing in green

Battery level 100%: Solid green

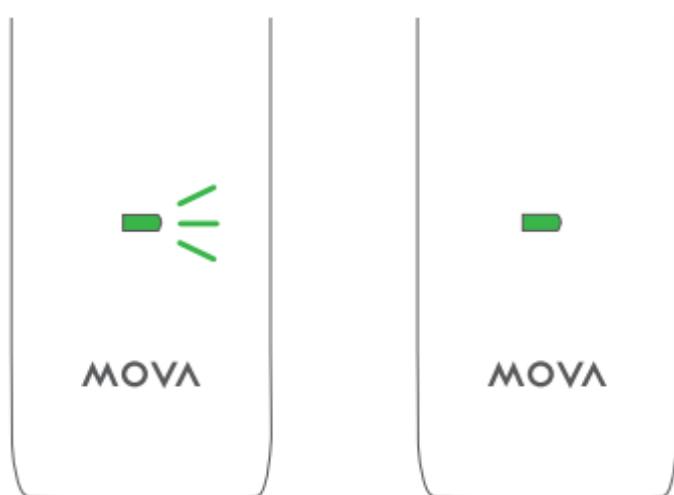Basic Parameters

Product name	Electric toothbrush		
Model	ATB13A		
Rated input voltage	5V	Rated input power	5W
Waterproof rating	IPX7	Battery capacity	900mAh
Standards	GB 4706.1-2005; GB4706.59-2008		

Names and Contents of Hazardous Substances

Part Name	Hazardous Substance					
	Lead (Pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent chromium (Cr(VI))	Polybrominated biphenyls (PBB)	Polybrominated diphenyl ethers (PBDE)
Brush head	○	○	○	○	○	○
Plastic parts	×	○	×	○	×	×
Hardware parts	×	×	×	○	×	×
Electronic parts	×	×	×	○	×	×
Power cord/accessories	×	×	×	○	×	×

This form is prepared in accordance with SJ/T11364.

○ It indicates that the content of the hazardous substance in all homogeneous materials of this component is below the limit specified in GB/T 26572.

×

It indicates that the content of the hazardous substance in at least one homogeneous material of this component exceeds the limit specified in GB/T 26572.

The number in this mark indicates that the product has an environmentally-responsible service life of 10 years under normal use.

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Manufacturer: Shaoxing Yueshi Technology Co., Ltd.