

SENCOR®

SOC 4110BL
SOC 4111GD



TOOTHBRUSH WITH INDUCTIVE CHARGING

Translation of the original manual




SENCOR®



SOC 4110BL
SOC 4111GD



A

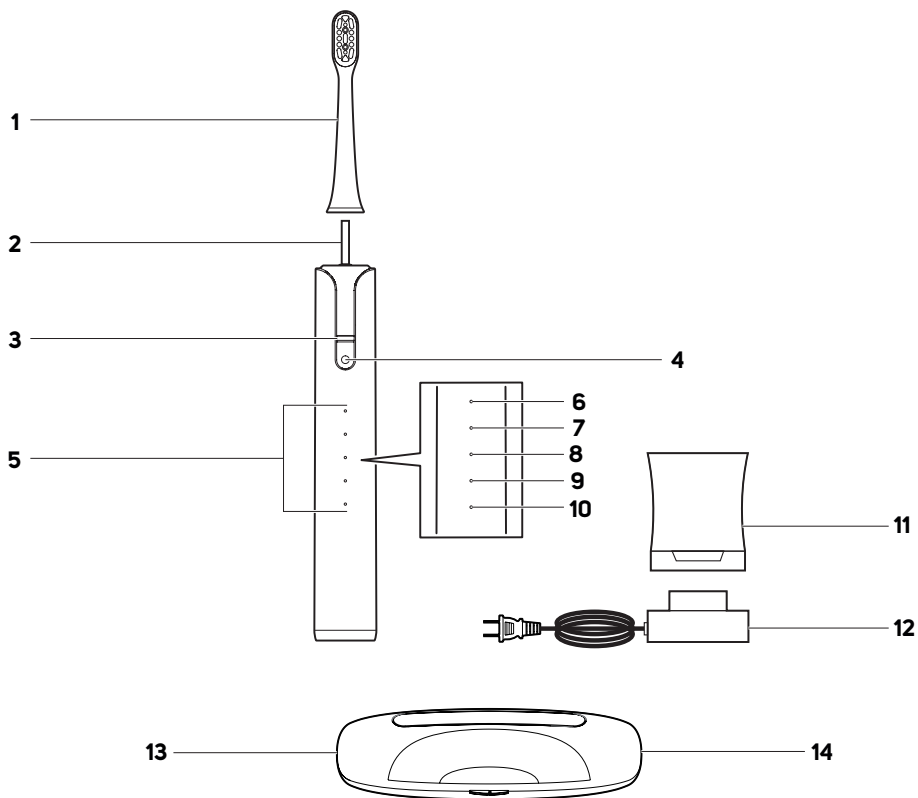

B

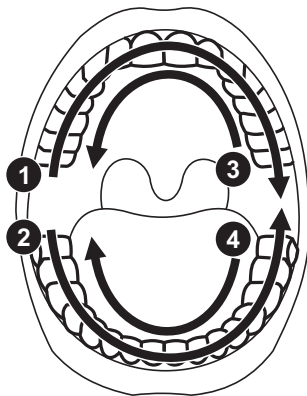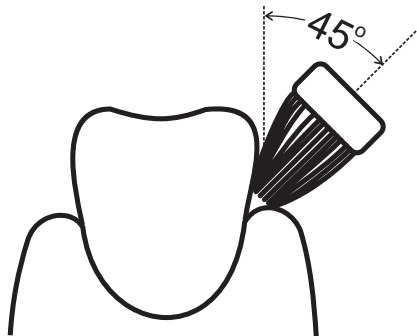

Important safety instructions

READ CAREFULLY AND STORE FOR FUTURE USE.

- Toothbrushes may be used by children and persons with reduced physical or mental capabilities or with a lack of experience and knowledge, if they are supervised or have been instructed to use the appliance in a safe manner and understand the potential dangers.
- Children must not play with the appliance. Cleaning and maintenance performed by the user must not be performed by unsupervised children.
- In the event that the power cord is damaged, have it repaired at a professional service centre to prevent causing a dangerous situation. It is forbidden to use the appliance if it has a damaged power cord.
- The parts of the appliance that need to be mounted or stood upright, must be mounted or stood upright in such a way that they cannot fall into water.
- Use the supplied original charging base to charge the toothbrush. Before connecting the charging base to a power socket, make sure the voltage stated on its rating label corresponds to the voltage in your power socket.

ELECTRICAL SAFETY

- Before connecting the charging base to an el. power socket, make sure the voltage stated on its rating label corresponds to the voltage in your power socket.
- Connect the appliance only to a properly grounded socket. Do not use an extension cord.
- Protect the base stand and its power cord against the ingress of water. Do not place it in the vicinity of vessels filled with water, e.g. sinks, baths, etc., do not spray it with water, submerge it in water or wash it under running water.
- Do not connect and disconnect the power cord of the charging base to and from the USB adapter with wet hands.
- Ensure that the USB power adapter plug into which the power cord of the base stand is connected does not come into contact with humidity, water or any other liquid. If precisely this happens, carefully dry the power plug before using it again.
- Do not disconnect the power cord from the toothbrush or from the USB power adapter by pulling on the power cord. This could damage the cord or the el. power socket. Disconnect the power cord from the power socket by pulling on the connectors.
- Do not place heavy items on the power cord. Make sure the power cord does not hang over the edge of a table and that it is not touching a hot surface or sharp objects.
- Keep the power cord dry.
- If the power cord or the USB power adapter is damaged, it must be replaced by the manufacturer, its service technician or by a similarly qualified person, as this will prevent a dangerous situation from arising. It is forbidden to use a damaged power cord or USB power adapter.
- Do not wind the power cord around the body of the toothbrush.
- The appliance is equipped with a built-in long-life battery that can only be removed at a professional service centre. Never attempt to remove it from the appliance yourself.



Attention:

Never disassemble the appliance for the purpose of removing the battery. Never short circuit the appliance or throw it into a fire. Incorrect use, risks releasing toxic substances or causing a fire or explosion.

- Do not repair the appliance yourself or make any adjustments to it. Have all repairs or adjustments of the appliance performed at an authorised service centre. Tampering with the appliance during the warranty period may void the warranty.
- To prevent the risk of injury by electrical shock, never submerge the power cord and the USB power adapter in water or any other liquid, or wash these parts under running water.
- If the charging base stand or the USB power adapter falls into water, under no circumstances continue using it, but rather hand it over to an authorised service centre for inspection.

SAFETY DURING USE

- This toothbrush serves to clean teeth and gums in the household. Do not use it for commercial purposes or for purposes other than those for which it is intended.
- Do not expose the toothbrush or its accessories to the effects of the outdoor environment, extreme temperatures, direct sunlight or locate it in the vicinity of an open flame, heating elements and other equipment or appliances that could be potential sources of heat.
- Charge, use and keep the toothbrush at room temperature.
- For hygiene reasons, toothbrush heads should be used by only a single user. Regularly replace the toothbrush head every 3 months, or sooner if it shows signs of wear. Use only original replacement toothbrush heads.
- Do not insert the toothbrush or its accessories into any bodily cavities, except for the mouth.
- When the toothbrush is used by children, supervise them while they are brushing their teeth until they acquire the proper tooth brushing technique.
- In the event that dental work or oral surgery was performed in the past two months, consult the use of a vibrating toothbrush with your dentist.
- In the event that gum bleeding occurs when the toothbrush is used and if bleeding does not stop even after a week of its use, contact your dentist.
- If your toothpaste contains peroxide, sodium bicarbonate or another bicarbonate (frequently used in whitening toothpastes), then thoroughly clean the brush and the handle after use.
- In the event that using the toothbrush causes unpleasant or painful feelings, stop using the toothbrush and consult this situation with your dentist.
- In the event that you have a pacemaker or another implanted device, consult the use of the toothbrush with your doctor or manufacturer of the implant.
- Perform cleaning and maintenance according to the instructions in chapter Cleaning and maintenance. Do not wash the toothbrush or its accessories in a dishwasher.
- Do not make any modifications to the toothbrush or its accessories. By tampering with the product you risk voiding your legal rights arising from unsatisfactory performance or quality warranty.

EXPLANATION OF SYMBOLS PLACES ON THE PRODUCT OR IN THE ACCOMPANYING DOCUMENTATION



This product meets all the basic requirements of the applicable EU directives.



This symbol on products or original documents means that used electric or electronic products must not be added to ordinary municipal waste.



DISPOSAL OF BATTERIES

Prior to disposing of the product or handing it over at a specialised collection facility, it is necessary to ensure that the battery is removed. Before removing the battery, it is necessary to start the vacuum cleaner and wait until it is completely flat. The battery contains environmentally damaging substances and therefore does not belong in standard household waste. Their correct disposal will be performed by the respective collection facility.



Safety class for electrical shock protection:

Class II – Electrical shock protection is provided by double or heavy-duty insulation.

.....

Toothbrush with inductive charging

User's manual

- Thank you for purchasing a SENCOR brand product and we hope that you will be happy with it.
- Prior to using this appliance, please read the user's manual thoroughly, even in cases, when one has already familiarised themselves with previous use of similar types of appliances. Only use the appliance in the manner described in this user's manual. Keep this user's manual in a safe place where it can be easily retrieved for future use. In the event that you hand this appliance over to somebody else, make sure to also include this user's manual.
- Carefully unpack the appliance and take care not to throw away any part of the packaging before you find all its parts. We recommend saving the original cardboard box, packaging material, purchase receipt and responsibility statement of the vendor or warranty card for at least the duration of the legal liability for unsatisfactory performance or quality. In the event of transportation, we recommend that you pack the appliance in the original box from the manufacturer.



Attention:

Do not place any metal items into the glass together with the toothbrush. They could affect the charging process and cause the charging base to overheat.



Note:

When charging the toothbrush for the first time after purchase or when it has not been used for more than three months, it is necessary to charge it for 24 hours.

The toothbrush may heat up during the charging process. This is normal and does not indicate a malfunction.

- When the charging process is complete, we recommend disconnecting the charging base from the power socket.
- The toothbrush is designed to regularly run for up to 3 weeks after being charged.
- Although it is possible to recharge the toothbrush after it is fully discharged, in order to maintain the optimal lifetime of the battery, we recommend charging it at regular intervals (e.g. after 2-3 weeks of use).

DESCRIPTION OF THE APPLIANCE

A1 Toothbrush head	A8 WHITE mode – whitening
A2 Connecting pin	A9 MASSAGE mode – gum massaging
A3 Charging indicator	A10 POLISH mode – for shiny teeth
Lit = battery is charging	A11 Glass
Not lit = battery is fully charged or is completely flat.	When placed on the charging base, it serves to charge the toothbrush
AĖ Main switch/cleaning mode toggle switch	A12 Charging base
A5 LED cleaning mode indicators	A13 USB-C slot for USB charging cable
A6 CLEAN mode – standard cleaning	A14 Travel case with UV filter
A7 SOFT mode – gentle cleaning mode for sensitive gums	

BEFORE FIRST USE

- Before first use, rinse the head of the toothbrush with clean warm water.
- Wipe the body of the toothbrush with a dry or a very lightly dampened wiping cloth. Do not submerge the charging base stand in water and do not wash it under running water.



Attention:

No part of this toothbrush or the charging base stand is intended for washing in a dishwasher!

- Prior to first use, charge the toothbrush battery according to the following chapter.

CHARGING USING THE CHARGING BASE

- This toothbrush is equipped with a contact contactless (inductive) charging system, which functions on the principle of energy transfer via an electromagnetic field. It is thus sufficient, to stand the toothbrush in the glass with the charging base and to connect the charging base to an el. power source.
- Place the glass with the charging base on an even and stable surface and connect the power cord to a power socket.
- Place the turned off toothbrush into the glass. Charging will commence automatically and the battery LED indicator on the body of the toothbrush will light up.
- When the battery is fully charged, the LED indicator will turn off.

CHARGING USING THE TRAVEL CASE

- The travel case protects the toothbrush and the accessories against loss, impact, scratching and soiling both when travelling and when at home. In addition to this, it is equipped with a USB-C slot, to which the supplied USB cable can be connected and the toothbrush battery can be recharged in this way if necessary.
- Place the toothbrush into the case.
- Plug the cable with the USB-C connector into the case. Connect the other end of the cable either to a compatible USB power adapter and then connect this to a power socket, or connect it to a USB port on a suitable device, such as a PC, notebook, etc.). Charging will commence automatically and the battery LED indicator on the body of the toothbrush will light up.
- When the battery is fully charged, the LED indicator will turn off. Disconnect the USB cable from the travel case and from the el. power source.

RECOMMENDATIONS FOR EXTENDING THE LIFETIME OF THE BATTERY

- In order for the battery to serve you as long as possible, charge, use and store the toothbrush in a dry environment at an ambient temperature of 15-35 °C. Lower or higher temperatures have a negative impact on the lifespan of the battery.
- When the toothbrush is not used for an extended period of time, charge it at least once every 6 months.

OPERATION

PREPARATION FOR BRUSHING

- Slide the toothbrush head on to the connecting pin on the toothbrush body. A small gap shall remain between the toothbrush unit and the ring on the bottom end of the head, which does not constitute a defect. To remove the toothbrush head, proceed in the reverse sequence.
- Rinse the head of the toothbrush with water and apply a small amount of toothpaste on to it.
- Place the toothbrush head into your mouth between the teeth and gum.
- Press the main switch to start the toothbrush.

SETTING THE BRUSHING MODE

- The toothbrush has five brushing modes:
CLEAN – mode for standard cleaning of healthy teeth;
SOFT – gentle mode for sensitive teeth or gums, we recommend using this mode when teeth are more sensitive, when tooth necks are exposed, periodontosis, implants and other dental substitutes,
WHITE – a mode suitable for whitening teeth
MASSAGE – pulse mode for gently massaging gums
POLISH – for shiny teeth
LED indicator of the selected mode will be lit.
- Toggle between these modes by repeatedly pressing the main switch:
-After first press of the main switch, the CLEAN mode will start.
-If you wish to select the subsequent mode, you must press the main switch again within 5 seconds of turning on the toothbrush. The LED indicator of the subsequent mode will be lit. If you press the main switch after more than 5 seconds of operation, the toothbrush will stop.



Note:
The toothbrush remembers the last used mode and therefore it is not necessary to reactivate it by pressing the mode toggle button. If you press the mode toggle button while cleaning your teeth, the toothbrush will switch to the new mode and the respective light will be lit. The toothbrush may heat up during the charging process. This is normal and does not indicate a malfunction.

CORRECT BRUSHING TECHNIQUE

- Place the toothbrush bristles on to the teeth at an angle of 45° in line with the gum, see figure B. Using short slow motions, move the toothbrush to the other teeth. Once the outer side of the upper right dental arch is brushed, move the toothbrush to the inner side of the dental arch. Again, proceed from back to front. Adjust the toothbrush angle to the tooth arrangement. When brushing the inner side of the front teeth, hold the toothbrush handle vertical so that the tip of the toothbrush head points upwards. Finally, brush the upper tooth surfaces.
- In the same way, brush the right half of the lower dental arch, the left half of the lower dental arch and finally also the left half of the upper dental arch. When brushing the inner side of the front teeth of the lower dental arch, hold the toothbrush handle vertical so that the tip of the toothbrush head points downwards.

BRUSHING TIMER

- This toothbrush has a timer function that momentarily pauses the toothbrush every 30 seconds until it turns off automatically, thereby informing the user that it is time to start brushing the next part of their teeth. This will ensure that all the dental quadrants are brushed for the same amount of time.

AUTOMATIC SHUT-OFF

- The toothbrush is equipped with a timer that will automatically stop it after 2 minutes of operation. This is the tooth brushing time recommended by dentists for the vibrating toothbrush.

ENDING THE BRUSHING PROCESS

- If you wish to stop the operation of the toothbrush before the automatic shut-off time limit, then repeatedly press the main switch until the toothbrush turns off.
- When you have finished using the toothbrush, clean it according to the instructions in chapter CLEANING AND MAINTENANCE.

CLEANING AND MAINTENANCE



Warning:
Prior to cleaning, the toothbrush must always be turned off and the charging base must be disconnected from the el. power source.
For cleaning the toothbrush and its parts, never use cleaning products with an abrasive effect, solvents, etc. Otherwise, its surface may be damaged. No part of this toothbrush or the charging base stand or its power cord is intended for washing in a dishwasher!

- For hygiene reasons, regularly clean the toothbrush after every use.
- Thoroughly rinse the toothbrush head under lukewarm running water to remove any toothpaste or other remains. If whitening pastes are used when brushing teeth, it is recommended to apply a drop of soap on to the toothbrush head and then to thoroughly wash it under running lukewarm water. Remove the toothbrush head from the handle at least once per week and clean it under running water in the location where it is attached to the handle.
- Clean the toothbrush handle using a lightly dampened cloth and then dry it.
- In the event that the charging base stand is soiled, wipe it with a clean dampened cloth and then wipe it dry. Ensure that the power cord plug does not come into contact with water or moisture.
- Never dry parts of the toothbrush over heat sources such as wood stoves, electric or gas stoves, etc.

REPLACING THE TOOTHBRUSH HEADS

- Regularly replace the toothbrush head every 3 months, or sooner if it shows signs of wear. Use only original replacement toothbrush heads. Replacement heads are available from you vendor both in online and brick&mortar stores.

STORAGE

- When not using the toothbrush for an extended period of time, place the washed toothbrush and brushing head in a dry, clean place that is out of reach of children. Keep the stored toothbrush and its parts away from extreme temperatures, direct sunlight and excessive humidity and do not store these in an overly dusty environment.

TECHNICAL SPECIFICATIONS

Toothbrush	
Battery.....	Li-Ion / 3.7 V / 800 mAh
Operating time when fully charged.....	up to 180 minutes
Noise level.....	<65 dB(A)
Protection class level against damaging entry of water.....	IPX7
Vibration speed.....	40,000 vibrations/minute
Charging base	
Input	100–240 V~, 50-60 Hz, 2.0 W
Protection class level against damaging entry of water.....	IPX7
Travel case	
Input.....	USB-C, 5 V , 0.6 A

The declared noise emission level of this appliance is lower than 65 dB(A), which represents a level A of sound power with respect to a reference sound power of 1 pW.

Protection class level against damaging entry of water:
IPX7 – The appliance is protected against damaging entry of water when submerged at a depth of 1 metre for 30 minutes.

INSTRUCTIONS AND INFORMATION REGARDING THE DISPOSAL OF USED PACKAGING MATERIALS

Dispose of used packaging material at a site designated for waste in your municipality.

DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT

For proper disposal, renewal and recycling, hand over these products at determined collection points. Alternatively, in some European Union states or other European countries you may return your products to the local retailer when buying an equivalent new product. Correct disposal of this product helps save valuable natural resources and prevents potential negative effects on the environment and human health, which could result from improper waste disposal. Ask your local authorities or collection facility for more details. In accordance with national regulations penalties may be imposed for the incorrect disposal of this type of waste.

For business entities in European Union states

If you want to dispose of electrical or electronic equipment, ask your retailer or supplier for the necessary information.

Disposal in other countries outside the European Union.

If you wish to dispose of this product, request the necessary information about the correct disposal method from the local council or from your retailer.

DISPOSAL OF USED BATTERIES

Prior to handing over the appliance at a specialised collection facility, it is necessary to ensure that the built-in battery is removed. Battery removal will be performed at the authorised service centre. In the event that the product has been disassembled for the purpose of ecological disposal then it must not be reassembled.

.....

Changes to text and technical parameters are reserved.

