The illustrations in this manual are for reference only.  
The actual appearance of the product may vary from what is shown.

# Table of Contents

EN	User Manual.....	001
PL	Instrukcja obsługi.....	009
SV	Användarhandbok.....	017
NO	Brukerhåndbok.....	025
FI	Käyttöohje.....	033
DK	Brugervejledning.....	041
RO	Manual de utilizare.....	049
AR	مدخستسما ليلد.....	057
HE	שמתשמל יודמ.....	065
ES	Manual de usuario.....	073
TR	Kullanım Kılavuzu.....	081
UA	Посібник користувача.....	089

## Safety Tips

### Safety Warnings

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and use maintenance shall not be made by children without supervision.
- The appliance is not intended for use by individuals (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure they do not play with the appliance.
- This appliance contains a non-replaceable battery. If the battery level remains consistently low, please seek professional assistance. For proper disposal, follow the battery removal instructions.
- Do not put the device or battery in fire or expose them to high temperatures to prevent battery leakage or overheating explosion.
- Please charge the product within a temperature range of 0°C to 45°C.
- Use an adapter with 5V  $\overline{=}$  output voltage to charge this product.
- Unauthorized user modification or repair of this product is strictly prohibited. If the product or power cord is damaged or battery leakage is detected, discontinue use immediately. To avoid danger, please let professionals of the manufacturer handle it.
- Do not use paint thinner, benzene, alcohol and other chemicals to clean this product.
- This product is a Class III appliance. The maximum measured output power of the power supply (power bank) is less than 15W. When connecting the power supply (power bank) to the charging cable for charging the appliance, its rated voltage should be 5V $\overline{=}$ , rated output current 1A, and it must simultaneously meet the requirements of Safety Extra-Low Voltage (SELV).

### Precautions for Electric Toothbrush

- When using this product for the first time, you may experience slight bleeding or tooth discomfort, as your gums have not yet adapted to the efficient cleaning. You can start with the Sensitive Mode or customize a gentler brushing mode. This phenomenon typically resolves within 1-2 weeks.
- If bleeding persists for more than 2 weeks, please stop using this product and consult your dentist.

## Safety Tips

### Precautions for Electric Toothbrush

- Users with severe periodontal disease or those undergoing dental repairs, orthodontics or implants should use this product under the guidance of a dentist.
- This product complies with safety standards for electromagnetic devices. For any questions related to pacemakers or other implanted devices, please consult your physician or the manufacturer of the implanted device before use.
- This product can only be used with the matching brush head.
- Using other brush heads may damage the tooth surface and gums, or affect the lifespan of the electric toothbrush. It is recommended to replace the brush head every 3 months to ensure the cleaning ability of the toothbrush and the cleanliness of the brush head.

### Maintenance

- Clean the brush head thoroughly after each use to prevent bacterial growth.
- Replace the brush head regularly; every 3 months is recommended.
- Avoid storing the toothbrush in a damp environment for extended periods to prevent bacteria buildup and battery damage.
- Avoid exposing the toothbrush to extreme heat or cold to protect battery life.
- Do not drop or impact the toothbrush to avoid internal damage.
- Regularly clean the casing and the connection area between the brush head and the body to prevent dirt buildup.
- If not in use for an extended period, clean all components and store the toothbrush in a cool, dry place away from direct sunlight.

### Battery Removal

When disposing of this product, the built-in rechargeable battery must be removed and recycled in accordance with local laws and regulations.

Dismantling the battery is not recommended under other circumstances.

1. Ensure the device is powered off and the battery is fully discharged before removing the battery.
2. Use a tool to open the bottom cover of the electric toothbrush and detach it from the main body.
3. Push inward on the latching points on both sides to remove the entire battery assembly.
4. Cut the connection cable between the battery and the motherboard body, take out the battery, and dispose of it at a government-designated location.

## Product Introduction

### Packing List

Electric toothbrush handle × 1, Whitening brush head × 1, Sensitive care brush head × 1, Toothbrush holder × 1, Type-C charging cable × 1, User manual × 1

### Part Description

#### Four-Zone Floating Brush Head

Each zone floats independently to adapt to the surface of your teeth. Shock absorption to protect gums without harming teeth

#### Four-Zone Ring Light

#### Power Button

#### Mode Light

#### Mode Button

#### Toothbrush Holder

#### Waterproof Cover

Keep it closed when not charging

#### Type-C Charging Port

#### Travel Case for Brush Head

Convenient brush head storage during travel


#### Type-C Cable

## Function and Usage

### Instructions for Use

- ① Insert the brush head into the handle.
- ② For first use, press the power button to exit flight mode.  
For daily use, pick up the toothbrush to activate.

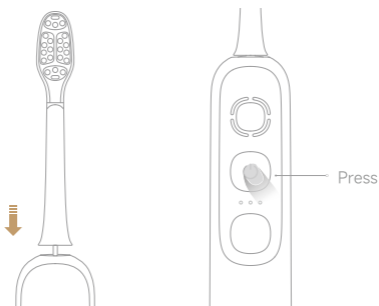- ③ **Start/Pause:** Short press the power button.  
**Switch mode:** Short press the mode button.  
\*In Exclusive Mode, you can customize the toothbrush's intensity with 10 levels. For specific instructions, please refer to the "Exclusive Mode" section.  
\*After 10 seconds of inactivity, the toothbrush will automatically enter sleep mode.

#### Sensitive Mode

Suitable for first-time users and those with sensitive teeth or gums

#### Cleaning Mode

Suitable for daily use and effective overall cleaning

#### Exclusive Mode

Suitable for individuals with specific requirements for toothbrush cleaning power

### Toothbrushing Coverage Reminder

- The four-zone light ring corresponds to the four tooth zones of your mouth: upper left, upper right, lower left and lower right. After brushing, the light bar corresponding to an insufficiently cleaned tooth area will turn red, while a fully cleaned tooth area will display white.  
\*Users can re-brush until all zones display white.

## Function and Usage

### Toothbrushing Coverage Reminder

Insufficient cleaning in the upper left zone

All areas cleaned

### Exclusive Mode

- ① Short press the Mode Button to switch to Exclusive Mode. Then press and hold the Mode Button to enter the custom adjustment state for Exclusive Mode.
- ② Short press the Mode Button to adjust the cleaning intensity from level 1 to 10.
- ③ After selecting the desired cleaning intensity level, press and hold the Mode Button or short press the Power Button. Once the mode indicator flashes three times, intensity customization of the Exclusive Mode is complete.

Long Press

Enter Exclusive Mode

Short Press

Adjust the cleaning intensity

Long Press

Complete intensity customization

### Flight Mode

- Press and hold the Power Button to enter Flight Mode.
- The toothbrush will automatically enter Flight Mode if it is shaken multiple times without being used.

## Function and Usage

### Flight Mode


Long Press  
Enter/Short Press Exit  
Flight Mode

Passively Enter  
Flight Mode

### Language Switch And Voice Control

- ① Hold the toothbrush with the brush head pointing downward at an angle of 60 to 120° relative to the horizontal plane.
- ② Press and hold the Mode Button to enter the setting state.
- ③ After each voice prompt, press and hold the Mode Button to switch to the next voice option.
- ④ During the voice prompt, short press any button to select the corresponding language and voice on/off setting.
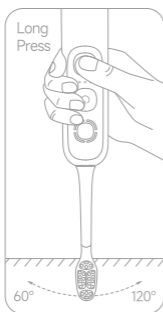
### Battery Level Alert

Mode Light Status	Toothbrush Battery Status	Working State
○ White light breathing	Charging	In progress


## Function and Usage

### Battery Level Alert

Mode Light Status	Toothbrush Battery Status	Working State
○ White light solid	Fully charged	Ready to work
● Red light solid	Low	Recharge soon
● Red light flashing	Critically low	Cannot operate—charge immediately

### Charging

- Open the waterproof cover at the bottom. If moisture is present at the charging port, wipe it dry.
- Gently open the waterproof cover at the bottom of the product. Do not pull it axially with force.
- Use the supporting Type-C charging cable and a power adapter or power bank that meets national standards. Do not use poor-quality chargers.
- Do not use the toothbrush while charging.
- Do not immerse the toothbrush in water or any other liquid while charging.
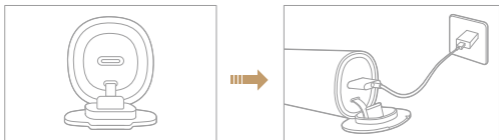
## Basic Parameters

Product name	MOVA Fresh 10 Electric Toothbrush		
Model	ATB22B		
Rated input voltage	5V $\overline{=}$	Rated input power	5 W
Waterproof rating	IPX8	Battery capacity	2000 mAh

We, XingKuang Innovation Technology (Suzhou) Co., Ltd., hereby declare that this equipment is in compliance in compliance with the applicable Directives and European Norms, and amendments. The full text of the EU declaration of conformity is available at the following internet address: [www.mova.tech](http://www.mova.tech)

This marking indicates that this product should not be disposed with other household wastes throughout the EU.

To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmentally safe recycling.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

This appliance contains batteries that are non-replaceable.

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Scan the QR code to  
download the MOVAhome APP  
Connect to MOVA electric  
toothbrush

Scan the QR code to view the  
product tutorial and electronic  
manual

User manual version number: V1.1

Service website: [www.mova.tech](http://www.mova.tech)

Manufacturer: XingKuang Innovation Technology (Suzhou) Co., Ltd.

Address: 12th Floor, Building 2, Wisdom Valley Park, Suzhou (Taihu)

Software Industrial Park, No. 1463 Wuzhong Avenue, Yuexi Street,

Wuzhong Economic Development Zone, Suzhou, Jiangsu Province,

P.R. China.

Made in China