3.2 Attaching the gel film

1. Carefully remove the blue protective film.

2. Carefully apply the gel film onto the EMS pad and then carefully peel away the protective film. Ensure that the edge of the gel film does not protrude over the EMS pad.

Gel film that has been applied slightly askew will not affect the function.

⚠ WARNING!

- Make sure that you attach the gel film to the device on the side where the blue protective film was, otherwise the skin could become irritated.

3.3 Usage

Stick the EMS pad onto the part of the body to be treated (see image).

1. Press and hold down the "+" button until you hear a long beep.

2. Press the "+" button again to start the Level 1 stimulation. Keep pressing the "+" button until you reach the desired intensity.

3. Press the "-" button to reduce the intensity. To turn off the EMS pad, press and hold down the "-" button until you hear a long beep. The device switches off automatically after 20 minutes.

Re-apply the protective film to the gel film after use.

For optimal results we recommend a 4-week training plan. Please note the following points:

1. 1-2 days rest.
2. Increase the intensity level on a daily basis.
3. A higher intensity level can lead to better results. However, the stimulation should never be uncomfortable.
4. Complete the full 20-minute program each time you use the device.

Wait at least 30 minutes after the stimulation programme is complete before using the device again.

3.4. Program

Cycle	3	2	3
Pulse width	400 µs	400 µs	400 µs
Pulse frequency	4 Hz	6 Hz	50 Hz
On-time	60 s	20 s	17 s/3"
Off-time	1.5 s	1.5 s	1.5/10 s**

On-time (s) = switch-on time/cycle in seconds (contraction). Off-time (s) = switch-off time/cycle in seconds (relaxation).

- * There are 3 individual rounds in cycle 3.
- ** The off-time is 10 seconds for every 3rd round in cycle 3.

Note: the product components do not require any maintenance.

5. Disposal

The empty, completely flat batteries should be disposed of through specially designated collection boxes, recycling points or electronic recycling. You are legally required to dispose of batteries. Note: the codes below are printed on batteries containing harmful substances: Pb = Battery contains lead, Cd = Battery contains cadmium, Hg = Battery contains mercury.

Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please contact your local authorities.

6. Technical specifications

Type EM 20 EMS Neck Pad
Output wave form Biphasic rectangular pulse
Pulse length 400 µs ± 10%
Pulse frequency 4-6-50 Hz ± 10%
Output voltage max. 70 Vpp ± 10% (1000 Ohm)
Output current max. 70 mApp ± 10% (1000 Ohm)
Voltage supply 3V CR2032 battery
Treatment time 20 minutes, the device then switches off automatically
Intensity Adjustable from 0 to 15


For optimal results we recommend a 4-week training plan. Please note the following points:

1. 1-2 days rest.
2. Increase the intensity level on a daily basis.
3. A higher intensity level can lead to better results. However, the stimulation should never be uncomfortable.
4. Complete the full 20-minute program each time you use the device.

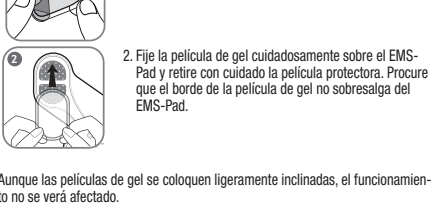Wait at least 30 minutes after the stimulation programme is complete before using the device again.

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